尺寸:
210x196mm

工艺要求:
70g双胶纸

半成品号:
E10769

MHK268 E10769

PUMP & COLLECT COMFORT COMBO - GEN. 1+ PUMP + PETAL COLLECTOR Product Manual English

Product Name: Pump & Collect Comfort Combo - Gen. 1+ Pump + Petal Collector
Materials: Pump - Food-Grade Silicone
Collector - Medical-Grade Silicone
Tube Cleaning Brush - Nylon, Stainless Steel
Dimensions: Pump - 84 x 57 x 147mm
Insert - 49 x 42 x 52mm
Collector - 110 x 110 x 41mm
Tube Cleaning Brush - 8 x 8 x 99mm
Pump Volume: 4oz. / 100ml

Product Description:

Meet Haakaa's Pump & Collect Comfort Combo, a carefully paired set designed for gentle pumping and discreet milk collection throughout the day. Whether actively expressing or catching light leakage while breastfeeding, this duo helps you collect milk comfortably wherever you are.

The Gen. 1+ Silicone Breast Pump uses natural suction for quiet, cord-free pumping, and now features the redesigned ComfortCARE™ Flange to help minimise friction on the breast. The included SpillGUARD™ Removable Insert helps keep milk safely contained during pumping for added confidence with every session!

The ultra-soft Petal Soft Silicone Milk Collector is designed to sit gently and discreetly inside your bra to catch let-down and light leakage without suction. Its slim, flexible shape naturally contours to the breast for comfortable everyday wear.

Made from soft, one-piece silicone designs that are simple to clean and can be sterilised for reuse again and again by simply boiling in water, this combo helps make milk collection feel easier, calmer and more convenient every day.

Features:

- Comfort-focused combo supports active pumping and passive milk collection for a relaxed breastfeeding experience.
- Express breastmilk with the pump using continuous natural suction – without cords or batteries!
- Redesigned ComfortCARE™ Flange features four ultra-soft ridges for a gentler, more secure fit on the breast.
- SpillGUARD™ Removable Insert helps reduce spills and leaks without interrupting pumping.
- Adjust suction by changing how much you squeeze the base, or fold back the flange for even stronger suction.
- The Petal discreetly collects let-down and light leakage without drawing out extra milk.
- Slim, petal-shaped collector is form-fitting for a subtle, natural look under everyday clothing.
- Internal silicone supports help the collector keep its shape and prevent leaks while moving.
- Pump & collector fit a wide range of breast shapes and sizes.
- Compact and lightweight for home, travel or everyday use.
- Pump's unique spillway & Petal's curved opening allow for quick milk transfer without spills.
- One-piece 100% silicone designs are safe in dishwashers, sterilisers and boiling water, with a dedicated cleaning brush included for thorough cleaning.
- BPA, PVC & phthalate-free.

⚠ WARNING:

- Check product condition regularly. Replace if this product shows any signs of damage.
- ALWAYS wash all parts of your breast pump and collector that come into contact with your breastmilk after every use.
- Cease use immediately and seek guidance from a lactation consultant if unable to express milk during two consecutive pumping session or if you experience pain or nipple damage.
- Never thaw frozen breastmilk in the microwave, directly in a pan on a stove element, or in boiling water. Do not use a microwave to warm stored breastmilk. Microwaving can destroy the nutrients in breastmilk and cause hotspots to develop that may cause seriousburns to your baby's mouth and skin.
- Using a breast pump or collector should not cause pain. If you have any red, hot, or inflamed areas on your breasts, discontinue use and consult a healthcare professional.
- If pain develops, or problems occur while using your breast pump, please consult your healthcare professional or lactation consultant.
- Do not use this pump during pregnancy unless approved by your maternal health professional, as pumping may induce labour.
- This is a single-user product. Use by more than one user without adequate cleaning may present a health risk and should

CAUTION:


- Read these instructions thoroughly to familiarise yourself with how the Haakaa Silicone Breast Pump and Petal Soft Silicone Milk Collector work before first use.
- You may find it helpful to practice on larger areas, such as your stomach, before your first time expressing milk.
- To avoid the risk of mould growth, ALWAYS thoroughly dry your breast pump and collector. Never store while wet.
- Do not use for any purpose other than expressing or collecting breastmilk.
- Pressing the flange too firmly to your breast may affect milk flow.
- Do not use suction strength any stronger than is comfortable, as this may damage breast tissue and affect your milk supply.
- Always ensure your nipple is properly centred within the flange during use.
- Your Haakaa Silicone Breast Pump may occasionally lose suction during use due to factors such as humidity, a full breast pump, or substances including moisturisers on your skin. Empty expressed milk, reposition the pump, or wipe excess skin products clean to rectify.
- This product and its packaging are not toys. Keep out of the reach of children and pets at all times.
- Adult use only.
- Do not use the Petal or pump while vigorously exercising or operating heavy machinery.
- Do not store near any sharp objects.
- Your silicone products may grow cloudy after repeated use, especially in hard water areas. This is normal and will not affect the performance of your product.
- As milk quantity increases in the pump, additional squeezing of the pump body may result in air or milk spurting out. Remove stopper, transfer expressed milk to a clean, sterilised container then reinsert stopper to continue pumping without leaks.
- The stopper is designed with small holes to stabilise pressure. To prevent leakage, avoid tipping pump upside down.

Cleaning and Care:


Wash before and after each use. Although dishwasher safe, we recommend hand washing with warm, soapy water and rinsing thoroughly. To clean the internal supports of Petal, use the included cleaning brush to scrub thoroughly with warm, soapy water and rinse well. To sterilise the pump, insert & collector, use a steam steriliser or boil in water for 2-3 minutes. Brush **NOT** suitable for microwaves, sterilisers or boiling water due to stainless steel components. Only soft bristle brushes or soft sponges should be used to clean this product as hard scourers may scratch the surface. Do not use any bleach-based cleaners or tablets to clean or sterilise your Haakaa products. Do not use UV sterilisers to sterilise this product. Store product in a cool, dry place and avoid direct sunlight.

Instructions:


Using the pump:

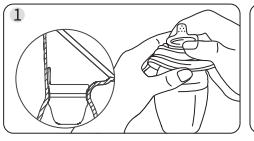


Wash your hands and all components before use. To insert the stopper, fold down the flange at the narrowest part of the neck. The protruding tab should face upwards and be aligned with the spillway. Fold flange back into place, ensuring there are no gaps.

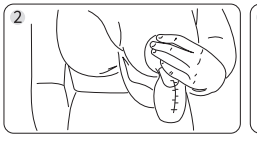


Place the pump flange over your breast, making sure your nipple is centred comfortably in the neck of the pump. You may need to adjust the position of your pump. Your nipple should not touch the sides of the pump or stopper.

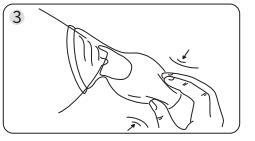

Once the pump is secure and comfortable, apply suction by squeezing the base of the pump to start the flow of milk.

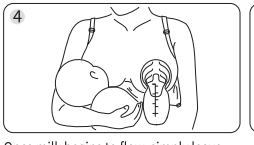

Once milk begins to flow, simply leave the pump alone to do its thing!

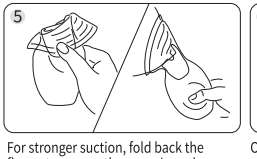

For stronger suction, fold back the flange to expose the pump's neck. Centre the opening over your nipple, squeeze the pump, attach to your breast as above, then return the flange to its original position.



Once you've finished expressing, remove stopper by lifting tab and pulling up before transferring your milk to a clean, sterilised container.


Using the Petal:



Hold the collector with the supports at the bottom, then centre the main opening over your nipple and hold to your breast.



Wear inside your bra to begin collecting let-down. The supports help minimise leaks during use so there is no need to press the back of your collector.

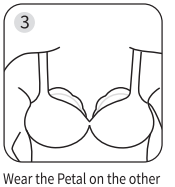

Wear the Petal on the other side while nursing to collect let-down or wear one on each breast to collect leakage.



The Petal can be worn while lying down – ensure Petal is secure inside your bra while lying down.

PLEASE NOTE: The Petal is suitable for collecting light daily milk leakage (holds up to 1 oz/30 ml). It will only collect milk that would otherwise have soaked into breast pads and will not draw out extra milk.

Empty collector regularly to prevent overflow and ensure your milk is stored safely.

TIPS: Encourage let-down by placing a warm flannel on top of your breast and massaging gently. While your child is feeding, use your Haakaa pump or Petal to the other breast to catch milk that would otherwise be lost in a nursing pad. Alternatively, use the pump to express milk, and wear the Petal on the other side to collect let-down. When the pump or Petal is full, empty into a clean, sterilised storage container or bottle and reposition.

The SpillGUARD™ Insert protects milk already in the pump body. If dropped, a few drops that haven't passed through the stopper may spill. This is normal.

NOTE: Always follow your country's guidelines and recommendations to safely express and store breastmilk. Volume measurements are an approximation and should be used for general reference only.

Every person is different and will produce different amounts of milk at different rates. Amounts may also vary depending on the time of day, or from one breast to the other.